
Two-way image splitter

The industry leading Optosplit II image splitter from Cairn Research divides an image into two separate, spatially identical components.


Displayed side by side on a single chip, this elegant device effectively doubles your research


capacity with one camera.



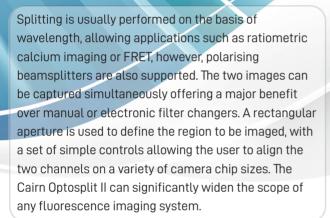


- Compact design with C-mount input and output ports as standard (F and T mount on request)
- Support for sensors up to 29.4mm diagonal
- Budget friendly alternative to two cameras
- 425nm to 875nm AR coatings on all optical surfaces
- 40mm diameter proprietary optics
- Simple & precise controls for image registration
- Interchangeable filter / dichroic holders
- Emission filter dimensions 25mm diameter
- Recommended dichroic dimensions 26x38x2mm (lambda/2 flatness)
- 1x, 1.3x and 1.7x magnification available
- Fixed or variable centre fully adjustable rectangular mask to delimit region of interest

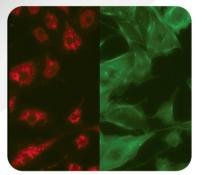

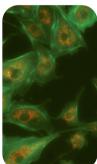

APPLICATIONS

- Förster Resonance Energy Transfer (FRET)
- Ratiometric calcium, voltage & pH imaging
- Simultaneous multi fluorescent probe imaging
- TIRF/Spinning disk confocal
- Polarisation studies (anisotropy)
- Simultaneous phase contrast / DIC and fluorescence
- Simultaneous dual Z depth imaging
- Single Plane Illumination Microscopy (SPIM)
- 3D super resolution PALM/STORM

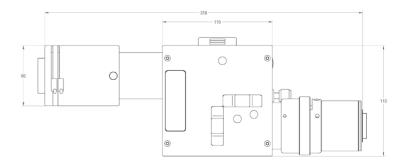



sales@cairn-research.co.uk +44(0)1795 590140 tech@cairn-research.co.uk

www.cairn-research.co.uk

MULTICHANNEL EMISSION SPLITTING RANGE

NO.1 IN OPTICAL PERFORMANCE, STABILITY AND USABILITY

DATASHEET



OptoSplit II & III

With an elegant configuration for simple side by side image splitting, the OptoSplit delivers high throughput imaging at a realistic price. Ideal for FRET, ratiometric imaging, polarisation studies and most simultaneous imaging applications requiring two or three images. User-configurable cubes and intuitive x, y adjustments offer convenience and simplicity.



Optosplit II Bypass

This builds on the success of the OptoSplit II, but adds a convenient single lever bypass mode making it more suitable for multi-user microscopes where simultaneous dual channel imaging is required for specific experiments alongside single wavelength recordings.



MultiSplit

Up to four channels simultaneously on one camera chip! The Multisplit uses the four quadrants of a single camera in a 2×2 square format. The Multisplit has the further possibility of simultaneous multi-depth imaging which is particularly attractive, as we can now do this at four depths rather than just two or three.



Multi Camera Adapters

Splitters for up to four channel imaging using multiple cameras (each). Perform simultaneous recording of multiple wavelengths, polarisation states or z depths using the full camera FOV without having to reduce their size. Variable rectangular aperture allows for the use of cropped sensor modes for the fastest speeds. Includes interchangeable camera mounts for C, F and T-mount cameras.



OptoMask

Enables precise FOV control for the high-speed, cropped sensor mode offered by several camera manufacturers including Andor and Teledyne..



OptoSpin 25 & 32

An intelligently designed, fast-spinning and stepping filter wheel, now with 25mm and 32mm wheels available and a range of microscope adapters. This slim unit has low inertia, enabling smooth operation and the ability to step between emission filters in 30ms, and spin continuously at 7500rpm when synchronised with a suitable light source. Change filters without moving the camera and mount two units together in the same 35mm optical path length for versatile combinations. (6 position for one filter wheel, 10 position for two).





